

1/1








P.I. Plug Cutter

D 544 420

4-1/2" with Obturator

3mm inside diameter, 5mm outside diameter

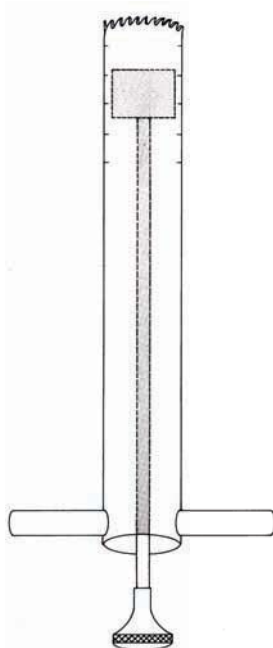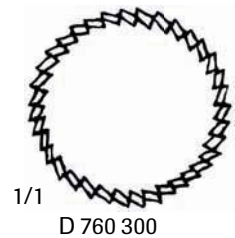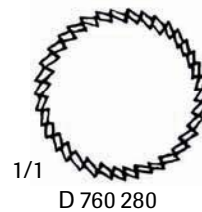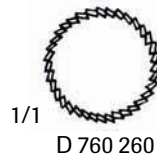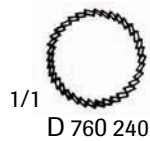

- These cutters are especially useful for arthrodesis of small and large joints by plug cutter graft technique.
- Plug cutters are designed to be used in pairs.
- Smaller diameter instrument will cut the plug opening and the larger one will harvest a bone graft that will fit into the opening.



P.I. Plug Cutter

6"

D 760 200 5mm inside diameter, 7mm outside diameter

D 760 220 7mm inside diameter, 9mm outside diameter

D 760 240 10mm inside diameter, 12mm outside diameter

D 760 260 12mm inside diameter, 14mm outside diameter

D 760 280 20mm inside diameter, 22mm outside diameter

D 760 300 22mm inside diameter, 24mm outside diameter